



**SRV947R  
REMOTE DISPLAY  
PLATFORM SCALE**

**Operating and Service  
Manual**

---

---

## TABLE OF CONTENTS

<b>TABLE OF FIGURES.....</b>	<b>2</b>
<b>PACKING CHECKLIST – MODEL SRV947R.....</b>	<b>3</b>
<b>ASSEMBLY.....</b>	<b>3</b>
<b>SYSTEM DESCRIPTION AND INTENDED USE.....</b>	<b>4</b>
<b>MAINTENANCE AND CLEANING.....</b>	<b>5</b>
<b>STORAGE AND TRANSPORTATION .....</b>	<b>5</b>
<b>SPECIFICATIONS.....</b>	<b>6</b>
<b>BUTTON FUNCTIONS .....</b>	<b>7</b>
<b>BASIC SYSTEM OPERATION.....</b>	<b>8</b>
<b>THEORY OF OPERATION .....</b>	<b>9</b>
<b>CALIBRATION.....</b>	<b>10</b>
<b>TROUBLESHOOTING .....</b>	<b>14</b>
<b>DISPLAY MOUNTING TEMPLATE.....</b>	<b>15</b>

---

---

## TABLE OF FIGURES

<b>Figure 1: SRV947R Button Label .....</b>	<b>7</b>
<b>Figure 2: Calibration Button .....</b>	<b>10</b>
<b>Figure 3: Select and Toggle Button .....</b>	<b>10</b>

---

---

## PACKING CHECKLIST – Model SRV947R Remote Display Platform Scale

√	DESCRIPTION	QUANTITY
	SCALE BASE	1 ea.
	SCALE DISPLAY	1 ea.
	MOUNTING HARDWARE	1 ea.
	AC WALL POWER SUPPLY (MEDICAL GRADE)	1 ea.
	CERTIFICATE OF CALIBRATION	1 ea.
	MANUAL	1 ea.

---

---

## ASSEMBLY

**STEP 1:** Unpack the scale system and check parts against the **PACKING CHECKLIST**. If there are any missing or damaged parts, call the Service Hotline at: 1-800-654-6360.

**STEP 2:** Determine installation location for scale and display, keeping in mind the display cable length is six (6) feet and the display power cord length is four (4) feet. The display may lay flat on a desk or end table. If desired, provisions for wall mount have been incorporated, and can be accomplished in the manner described below.

- A. The display should be installed by qualified personnel. The mounting hardware provided will cover most wall applications.
- B. Determine a position on the wall comfortable to the user. Mark the location of the two (2) mounting holes on the wall (holes are spaced 2.5” on center apart) using the template provided on page 15.
- C. Drill the two (2) pilot holes, no larger than ¼” diameter in the locations marked. Install the two (2) self-drilling screw anchors.
- D. Install the two (2) #8 x 1” flat head screws supplied, leaving approximately ¼” projecting from the wall.

**STEP 3:** Connect signal cable to the underside of the platform.

**STEP 4:** Place SRV947R scale system on a flat, hard floor surface. Ensure all six feet of the scale are in contact with the floor and scale does not rock.

**NOTE:** Scale should not be placed on carpet or uneven floor surfaces as this may interfere with scale accuracy.

**STEP 5:** Connect the display cable to the scale, then connect the power cable to the display and plug power adapter into wall.

---

---

## SYSTEM DESCRIPTION AND INTENDED USE


### SYSTEM DESCRIPTION

The SRV947R Portable Platform Scale is a low profile design, 2 in. height, which allows the animal easy step access. The SRV947R Scale can be used in a fixed location or as a portable scale. The Scale employs the latest in microprocessor and load cell technology to provide accurate and repeatable weight data. Six (6) identically matched transducers are strategically placed to ensure an accurate representation of the patient's weight.

The patient's weight is displayed on an LCD screen. With a push of a button, weight data may be viewed in either pounds or kilograms.

### INTENDED USE

The SRV947R Portable Platform Scale is a lightweight portable scale designed for use as a weighing system for animals not exceeding 400 pounds or 182 kilograms, gross weight. Please read the following warnings to avoid injury to the animal or attendant.

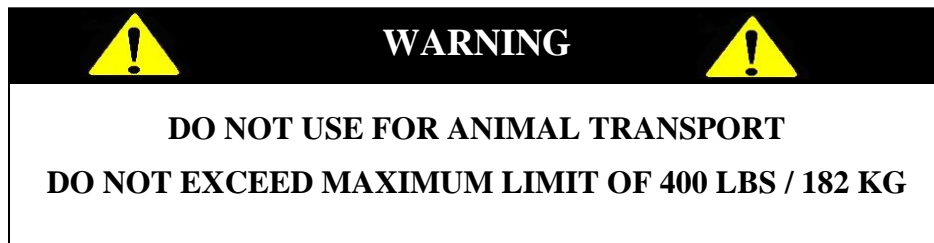


---

---

## MAINTENANCE and CLEANING

The SRV947R Portable Platform Scale is made with a polyurea-coated aluminum platform. Exercise caution when cleaning the display window as it is made of a clear plastic and can be scratched by abrasive cleaners. Mild soap and water is recommended for general cleaning and disinfecting.

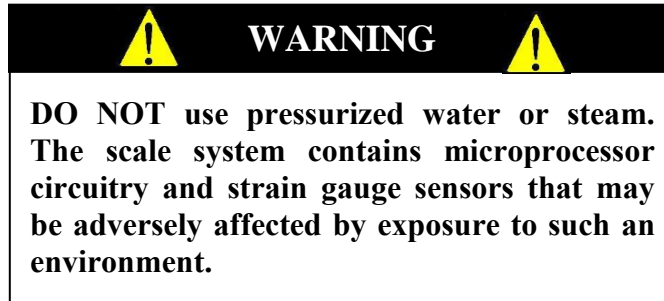

---

---

## STORAGE and TRANSPORTATION

### STORAGE

If storing this equipment for periods longer than three (3) months, remove the batteries. To maintain proper operation of this instrument, storage and transport conditions should not vary outside the following conditions: Relative Humidity 0% to 85%, Ambient Temperature 14°F to 122°F (-10°C to +50°C).

### TRANSPORTATION

The SRV947R Portable Platform Scale has four (4) recessed handles on the underside of the scale for easy portability.

## SPECIFICATIONS

<b>MAXIMUM WEIGHT CAPACITY</b>	400 lb / 182 kg
<b>PLATFORM SIZE</b>	24 in x 48 in x 2 in (61 cm x 122 cm x 5 cm)
<b>DISPLAY TYPE</b>	1.0" LCD display with 5 active digits
<b>DISPLAY RESOLUTION</b>	0.2 lb / 0.1 kg
<b>ACCURACY</b>	0.2% +/- 2 digits of displayed resolution for calibrated range
<b>AUTO ZERO</b>	One button operation
<b>AUTO WEIGH</b>	Scale will automatically turn on and lock in on a stable weight.
<b>AVERAGING</b>	Automatic digital filter
<b>POWER SUPPLY</b>	AC Wall Power Supply (Medical Grade)
<b>CALIBRATION</b>	Calibration is traceable to NIST standards
<b>OPERATING CONDITIONS</b>	Normal operating conditions for this product: Ambient Temperature Range: 68°F to 85°F (20°C to 30°C), Relative Humidity Range: 0%-85%. Avoid exposure to high-pressure water or steam.
<b>TRANSPORT and STORAGE</b>	Storage and transport conditions should not vary outside the following conditions: Relative Humidity 0% to 85%, Ambient Temperature 14°F to 122°F (-10°C to +50°C). Remove batteries if storing longer than three (3) months.

---

## BUTTON FUNCTIONS

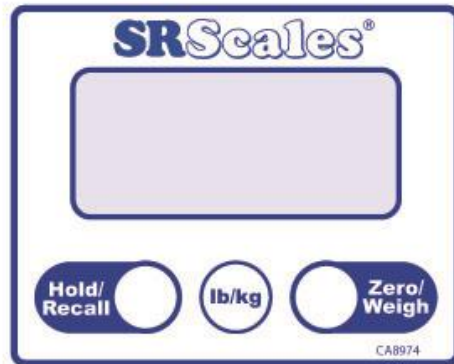


Figure 1: SRV947R Button Label

- ZERO** The “ZERO/WEIGH” button is used to clear the Recall memory and start a new weighing cycle. Press and hold for 3 seconds to zero the display.
- MODE** The “MODE” button is used to select pounds or kilograms.
- HOLD/  
RECAL** The “HOLD/RECALL” button is used to “lock in” the weight currently displayed, or to recall the last saved weight, if nothing is on the scale platform.

---

---

## BASIC SYSTEM OPERATION

### TURNING SYSTEM ON

Make sure the scale is free and clear of any obstructions and all six feet are touching a level surface. Press the “**ZERO**” button and hold for 3 seconds. The display will read “**ZERO**” then display “**0.0**” when system is ready for weighing.

### UNIT

Press the “**MODE**” button to toggle between pounds or kilograms. The system will remain in the same mode as last used unless manually changed.

### WEIGH

The “**WEIGH**” button is used to clear the “Recall” memory and start a new weighing cycle.

### HOLD/RECALL

The “**HOLD/RECALL**” button is used to “lock in” and store the weight currently displayed in system memory. With no weight on the platform pressing the button again will recall the last stored weight and display it for 15 seconds.

### AUTOHOLD

With autohold enabled, when the system detects a stable weight, that value will be stored in system memory. This value can later be recalled by pressing the “**HOLD/RECALL**” button with no weight on the platform.

---

### SCREEN DISPLAYED SYMBOLS



Scale is zeroed and stabilized

**MEM**

Weight is locked in memory (next stable weight will override memory).



Applied weight load has stabilized



---

---

## THEORY OF OPERATION

SR Instruments patient weighing systems are digital scales. Strain-gauge force cells convert the force of an applied weight into an analog signal. This signal is amplified by an operational amplifier and converted to a digital signal by analog to digital converter. The digital signal is transferred to a micro-controller where it is filtered, converted to appropriate units, and displayed on a liquid crystal display.

Strain-gauge force cells each contain four strain gauges mounted in a full Whetstone-bridge configuration. These bridges convert the physical movement of the force cell, due to the applied mass on the system, into minute changes in electrical resistance. These changes in resistance produce a voltage difference across the Whetstone-bridge, which is amplified by the operational amplifier.

The output of the operational amplifier is digitized by the analog to digital converter. The converter integrates the analog signal onto the integrating capacitor over a short interval. The integrating capacitor is then discharged at a rate proportional to the reference voltage applied to the converter. The residual voltage on the integrating capacitor is then multiplied by a factor and again discharged at a rate proportional to the reference voltage. The residual voltage from this discharge is again multiplied by a factor and again discharged. The time taken to discharge the capacitor is proportional to the voltage from the operational amplifier, which is proportional to the applied load on the force cells. The time is stored as a binary number in the analog to digital converter and is transferred to the micro-controller when the conversion is complete.

The micro-controller averages and filters the digital output of the analog to digital converter, subtracts the value saved during the system zero operation and scales the filtered output, then displays the result on the liquid crystal display. The micro-controller performs a rolling average of data for continuous weigh and averages the data before locking in on the reading.

## CALIBRATION



### IMPORTANT



**CALIBRATION CHECK** Qualified service personnel only should perform this procedure. Load cells have no user serviceable components and should not be tampered with for any reason. Re-calibration is generally not required, but should be verified periodically to ensure accuracy. The recommendation for calibration check is at least once every 12 months, or as individual maintenance policy requires.



### CAUTION



The integrated circuits and semiconductors on the printed circuit boards may be damaged by electrostatic discharge (ESD). Be sure to use proper handling precautions at all times.

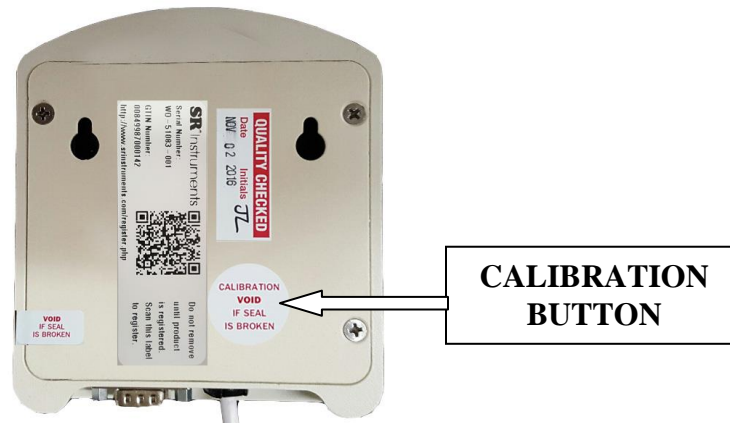


Figure 2: Calibration Button

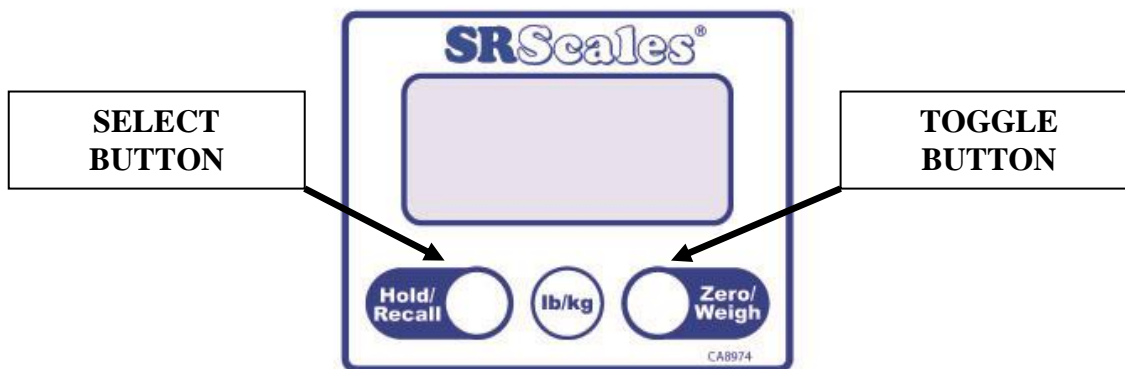


Figure 3: Select and Toggle Button

---

## INITIAL SYSTEM SETUP

When initially set up, calibration is factory set and re-calibration is not needed. The Local Gravitational Acceleration may have to be re-set for the current geographical location for improved accuracy. Procedure is found below.

## ACCESS SYSTEM SETUP

To access all settings below, break the Calibration Seal. Press the calibration button (Figure 2) using the pointed-end of a pen or similar item. CAL mode is now active.

Refer to Figure 3. The Right Button allows you to toggle through the menus without saving any parameters. The Left Button permits a menu item to be selected or a parameter saved.

**NOTE:** Ensure that nothing is in contact with the scale system during this procedure. Remove hands from the system when noting the displayed calibration results.

## CAL MENU

- “**data**” Displays saved calibration and model information.
- “**unit**” Set the available units; lb only, kg only, or selection of lb or kg with mode button.
- “**sec**” Set Auto Power Down Timer between 15 and 900 seconds. Setting is non-functional on this model
- “**grav**” Set Gravitational Acceleration in  $m/s^2$
- “**Full**” Applied calibrated mass at full scale.
- “**Half**” Used with three point calibration curve. Applied calibrated mass at a second point between zero and full scale, usually at half.
- “**zero**” No calibrated weight applied.
- “**2 pt**” Two-point calibration curve, zero and full scale.
- “**3 pt**” Three-point calibration curve; zero, half and full scale.
- “**save**” Saves the calibration data.
- “**quit**” Quits without saving the calibration data.

## Enter CAL Mode

To enter the CAL mode, press the calibration button (Figure 2) using the pointed-end of a pen or similar item. Toggle through the menu with Right Button (Figure 3). Select or Save data with the Left Button (Figure 3).

---

### Set the Gravitational Acceleration (Optional)

**NOTE:** Setting the local gravitational acceleration, “**gLOC**”, will not affect the calibration. The local gravitational acceleration, “**gLOC**”, will become the calibrated gravitational acceleration, “**gCAL**” only when it is recalibrated.

- a) Select the Gravitational Acceleration from the CAL menu by pressing the Left Button when the display shows “**grav**”.
- b) The display will show only the four digits to the right of the decimal point <9.XXXX>.
- c) Set the local gravitational acceleration in  $m/s^2$  by selecting the first digit and toggle through the digits using the Right Button until desired value is shown. Press the Left Button and continue the same with the remaining digits.
- d) Press the Left Button to save or Right Button to toggle to “quit”. Select “quit” by pressing the Left Button to quit without saving changes.
- e) Press the Left Button again to save and the display will show “**SAV’d**”.

### Set the Full-Scale value

- a) Select the Full-Scale from the CAL menu by pressing the Left Button when the display shows “**FULL**”.
- b) The display will show only the five digits <XXXX.X> with the thousand-place flashing <1>. Select past to the next digit by pressing the Left Button.
- c) Once in the hundreds-place, set the full scale value (maximum load) by toggling using the Right Button through the digits until desired value is shown. Select the desired value by pressing the Left Button and continue the same with the remaining digits.
- d) Press the Left Button to save or Right Button to toggle to “quit”. Select “quit” by pressing the Left Button to quit without saving changes.
- e) Press the Left Button again to save and the display will show “**SAV’d**”.

---

### Set the Half-Scale value

**NOTE:** The Half-Scale value is used only with the three-point calibration. It is a value between zero and the Full Scale values. It is usually close to half of the Full Scale value.

- a) Select the Half-Scale from the CAL menu by pressing the Left Button when the display shows “**HALF**”.
- b) The display will show only the five digits <XXXX.X> with the thousand-place flashing <1>. Select past to the next digit by pressing the Left Button.
- c) Once in the hundreds-place, set the half scale value by toggling using the Right Button through the digits until desired value is shown. Select the desired value by pressing the Left Button and continue the same with the remaining digits.
- d) Press the Left Button to save or Right Button to toggle to “quit”. Select “quit” by pressing the Left Button to quit without saving changes.
- e) Press the Left Button again to save and the display will show “**SAV’d**”.

### Calibrate using the Three-Point Calibration Curve; zero, half scale, and full scale.

**NOTE:** The scale can also be calibrated with a Two-point slope “**2Pt**”, although it will not be as accurate.

- a) Select Three-point calibration curve by pressing the Left Button when the display shows “**3Pt**”
- b) Apply no load and ensure stability. Zero the scale by pressing the Left Button when the display shows “**zero**”.
- c) Display will show “**HALF**”, apply the predetermined half scale value on the scale and press the Left Button.
- d) Display will show “**FULL**”, apply the predetermined full scale value on the scale and press the Left Button.
- e) Display will show “**SAVE**”, press the Left Button to save, display will show “**SAV’d**” or Right Button to toggle to “quit”. Select “quit” by pressing the Left Button to quit without saving changes.
- f) Save the parameter changes while in the CAL mode by pressing the Left Button when the display shows “**SAVE**”. To exit the CAL mode without saving the changes press the Left Button when the display shows “**quit**”.
- g) Press the calibration button (Figure 2) using the pointed-end of a pen or similar item.
- h) Place a new calibration seal over the calibration button access hole.

**NOTE:** Calibrated weight displayed should be within tolerance as per Model being calibrated (see “**SPECIFICATIONS**” page 6).

---

---

## TROUBLESHOOTING


SYMPTOM	REASON/CORRECTIVE ACTION
Weight reading is much lower than expected.	Check that the platform is clean underneath, stand on each of the four corners to see if one corner is not weighing correctly.
System fails to perform correctly	Check power supply connections. Make sure patient is clear of any obstacles and that nothing, including the caregiver, is touching the patient or scale. Set scale on smooth, hard floor surface.
<b>For additional information or assistance, telephone our Service Hotline: 1-800-654-6360 or e-mail: <a href="mailto:sri@srinstruments.com">sri@srinstruments.com</a></b>	

---

## DISPLAY MOUNTING TEMPLATE

Determine installation location for scale and display, keeping in mind the display cable length is six (6) feet and the display power cord length is four (4) feet. The display may lay flat on a desk or end table. If desired, provisions for wall mount have been incorporated, and can be accomplished in the manner described below.

- A. The display should be installed by qualified personnel. The mounting hardware provided will cover most wall applications.
- B. Determine a position on the wall comfortable to the user. Mark the location of the two (2) mounting holes on the wall (holes are spaced 2.5" on center apart) using the template provide below.
- C. Drill the two (2) pilot holes, no larger than 1/4" diameter in the locations marked. Install the two (2) self-drilling screw anchors.
- D. Install the two (2) #8 x 1" flat head screws supplied, leaving approximately 1/4" projecting from the wall.

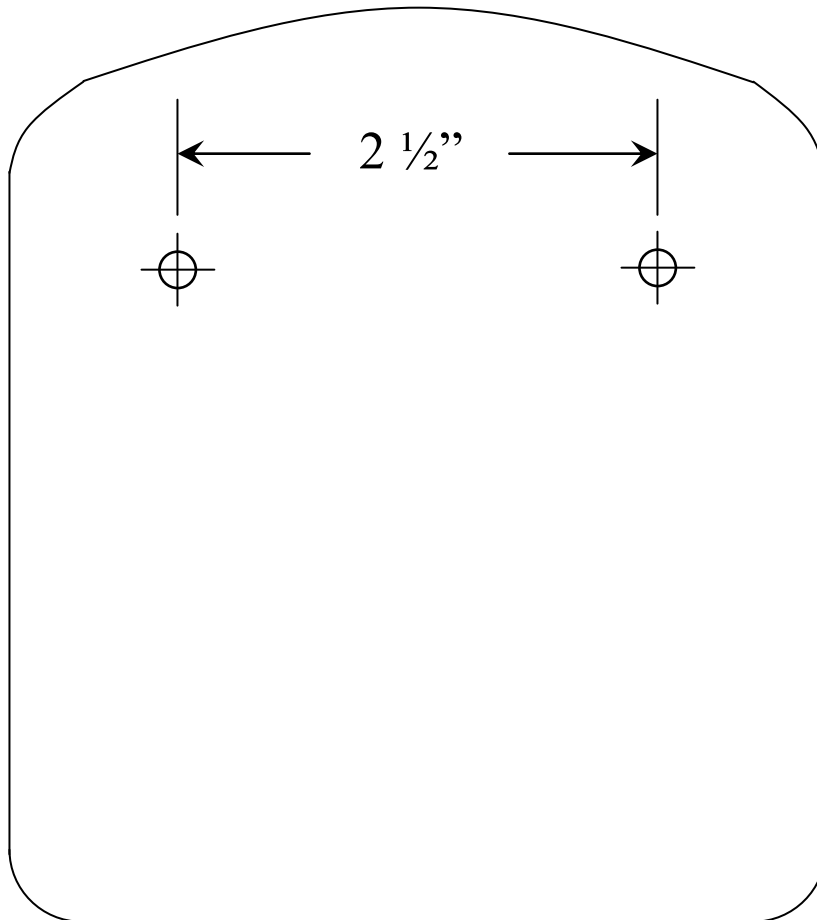



**Precision & Technology in  
Perfect Balance™**