
SPECIFICATIONS

Max Capacity:	200g
Min Capacity:	0.5g
Division:	0.1g
Power source:	4-“AA” batteries
Accuracy:	0.2% +/- 1 digit of displayed resolution for calibrated range

FEATURES

Portability:	Small and compact
Tare:	Single button
Auto power off:	120 seconds
Warranty:	Two (2) years

BATTERY INSERTION/ REPLACEMENT

1. Turn scale over. Remove the Battery Compartment Cover by pressing the tab towards the cover and lift it off.
2. Place four (4) – 1.5-Volt batteries of AA size into the battery compartment with correct polarity.
3. Replace the Battery Compartment Cover.

TO OPERATE SCALE

1. Place scale on a flat surface in a sheltered area.
2. Press the “ON” button to turn the scale on.

USING THE TARE FUNCTION

1. Place the scale on a flat surface in a sheltered area.
2. Press the “ON” button to turn the scale on.
3. When “0” appears on the display, place the empty container onto the weighing platform.

4. Press the “TARE” button and the display will return to zero, eliminating the weight of the container.
5. A “o” will appear in the lower left corner of the display, which indicates that the TARE feature is functioning.
6. Place the object to be weighed into the container and obtain its net weight.

MANUAL POWER OFF

After use, turn the scale off by pressing the “OFF” button to prolong battery life

AUTO POWER OFF FEATURE

The scale is also equipped with an “AUTO POWER OFF” function. The scale will turn off in approximately 120 seconds.

CAUTION TO BE EXERCISED

1. Place the scale in a sheltered area when in use. Air movement on the scale may affect the accuracy.
2. Always place object onto the platform gently. **DO NOT** weigh objects exceeding the full capacity as this may permanently damage the scale. . An “O-Ld” message will appear on the display if the scale is over-loaded.
3. Replace batteries when the battery symbol appears in the display when scale is turn on. When batteries are extremely low, “LO” will appear and unit will shut off.

CALIBRATION METHOD

1. With the scale **OFF**, **PRESS** the “**ON**” button five (5) times and hold. Hold until a random number appears.
2. Press the “**ON**” button. The display will read “**CAL**”.
3. After 2 seconds, the display will read 200 grams (do not use any other size weight as the calibration will be altered).
4. Place the 200 grams weight on the scale and press the “**ON**” button. The display will read “**PASS**”. Calibration is complete.

TWO (2) YEAR LIMITED WARRANTY

Each **SR Scales**® system is manufactured with high quality components. SR Instruments, Inc. warrants that all new equipment will be free from defects in material or workmanship, under normal use and service, for a period of two (2) years from the date of purchase by the original purchaser. Normal wear and tear, injury by natural forces, user neglect, and purposeful destruction are not covered by this warranty. Warranty service must be performed by the factory or an authorized repair station. Service provided on equipment returned to the factory or authorized repair station includes labor to replace defective parts. Goods returned must be shipped with transportation and/or broker charges prepaid. SR Instruments, Inc.’s obligation is limited to replacement of parts that have been so returned and are disclosed to SR Instruments, Inc.’s satisfaction to be defective. The provisions of this warranty clause are in lieu of all other warranties, expressed or implied, and of all other obligations or liabilities on SR Instruments, Inc.’s part, and it neither assumes nor authorizes any other person to assume for SR Instruments, Inc. any other liabilities in connection with the sale of said articles. In no event shall SR Instruments, Inc. be liable for any subsequent or special damages. Any misuse, improper installation, or tampering, shall void this warranty.

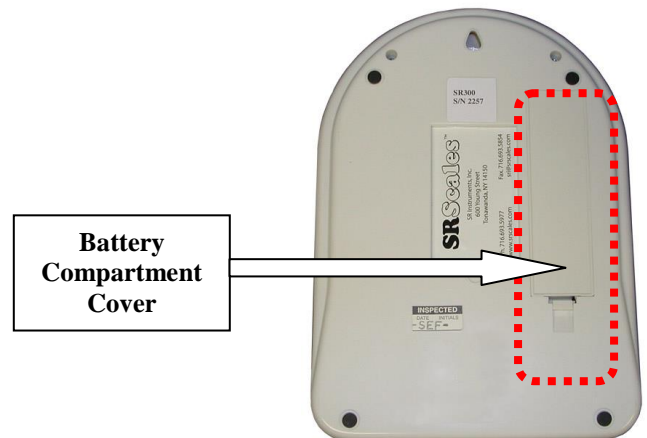
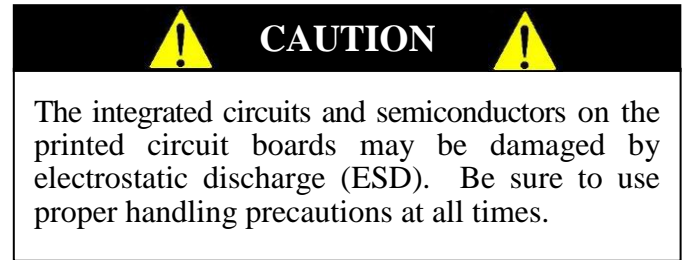


Figure 1: Battery Compartment Cover