



SR 735i

Wheelchair Scale

Operating and Service Manual

TABLE OF CONTENTS

TABLE OF FIGURES.....	2
PACKING CHECKLIST.....	3
ASSEMBLY.....	4
REPLACEMENT PARTS AND ACCESSORIES	5
SYSTEM DESCRIPTION AND INTENDED USE.....	6
MAINTENANCE AND CLEANING.....	6
STORAGE.....	7
SPECIFICATIONS.....	7
BUTTON FUNCTIONS	8
BASIC SYSTEM OPERATION.....	10
BASIC SYSTEM OPERATION CON'T	11
BATTERY REPLACEMENT	11
THEORY OF OPERATION	12
CALIBRATION.....	13
CONFIGURING SEND BUTTON.....	15
TROUBLESHOOTING	16
WARRANTY.....	17

TABLE OF FIGURES

Figure 1: Attach Wall Mounting Plate.....	4
Figure 2: Wall Mount Hole Locations	4
Figure 3: Install “D” Cell Batteries	4
Figure 4: Printer Mounting Bracket	5
Figure 5: Printer Display Connector.....	5
Figure 6: Button Display	8
Figure 7: Battery Replacement.....	11
Figure 8: Calibration Tolerance Table	13
Figure 9: Calibration Switch Diagram	13

PACKING CHECKLIST

SR735i Platform Scale

√	DESCRIPTION	QUANTITY
	Box 1	
	WHEELCHAIR PLATFORM BASE	1 ea
	DISPLAY UNIT	1 ea
	PACKAGE OF SIX (6) "D" CELL BATTERIES	1 ea
	CALIBRATION CERTIFICATE	1 ea
	WARRANTY CARD	1 ea
	MANUAL	1 ea

ASSEMBLY

STEP 1: Unpack the scale system and check parts against the **PACKING CHECKLIST**. If there are any missing or damaged parts, please call the service hotline at 1-800-654-6360.

STEP 2: Verify that the serial number on the Display Unit matches that on the Wheelchair Scale Assembly.

STEP 3: The SR735i Display may be wall mounted (**Figure 1**) using the Wall Mounting Plate. **Note:** See **Figure 2** for locating holes for wall mount application.

Figure 1: Attach Wall Mounting Plate

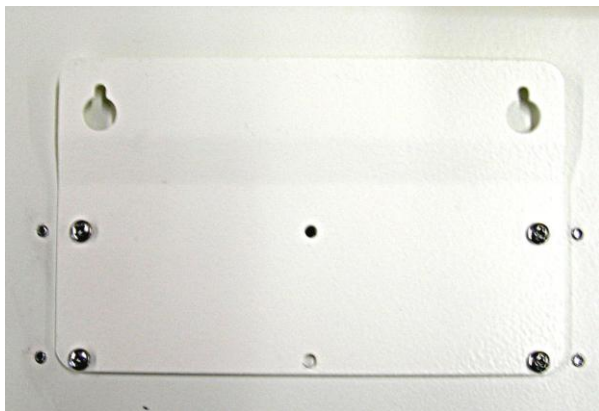


Figure 2: Wall Mount Hole Locations

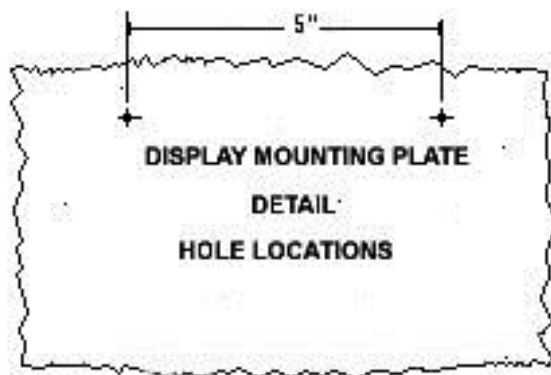


Figure 3: Install “D” Cell Batteries



STEP 4: (**Figure 3**) Remove the battery compartment cover by unscrewing the battery cover thumb screw and install the six (6) “D” cell batteries (**Figure 3**). Replace cover and tighten screw.

Continued next page

ASSEMBLY cont'd

INSTALLING OPTIONAL PRINTER

STEP 1: Mount the optional printer within distance of the display to allow connection of the printer cable. Use the three (3) holes in the Printer Mounting Bracket as a template to locate the holes needed to mount the Printer Mounting Bracket (**Figure 4**). Use the wall mounting hardware provided for installing the bracket.

Figure 4: Printer Mounting Bracket



Figure 5: Printer Display Connector



STEP 2: Carefully attach the printer to the Printer Mounting Bracket. Connect the printer cable from the printer to the connector on the bottom of the display (**Figure 5**). Tighten the two (2) printer cable screws to secure the cable to the connector.

REPLACEMENT PARTS and ACCESSORIES

Part #	Description
SR8327	Printer Kit (includes all printer accessories listed below)
FRAP1300	Printer
FRAP1300BR-01	Printer Mounting Bracket
FRAP1300-D9	Printer Cable
FRAP1300PCC	AAA Battery Holder
FRAP1300BP01	Rechargeable Battery Pack
FRBC1300	Printer Power Cable
FRTPT130012C	Paper, thermal printer 58mm (10 rolls)

SYSTEM DESCRIPTION and INTENDED USE

SYSTEM DESCRIPTION

The SR735i Wheelchair Scale System employs the latest in microprocessor and load cell technology to provide accurate and repeatable weight data. Four (4) identically matched transducers are strategically placed to ensure an accurate representation of the patient's weight.

The SR735i low power microprocessor circuitry is powered six (6) common "D" cell batteries, which will provide up to 400 hours of weight readings before needing replacement.

The patient's weight is displayed on a 16-character dot matrix LCD. The weight data may be viewed either in pounds or kilograms with a displayed resolution of 0.1 for each.




INTENDED USE

The SR735i Wheelchair Scale System is specifically designed for weighing both motorized-wheelchair and caregiver-assisted wheelchair bound patients and is a preferred means of gathering patient weight data of up to 1000 pounds or 454 kilograms.



MAINTENANCE and CLEANING

SR735i Platform Scale Display is made of powder-coated aluminum. Exercise caution when cleaning the display window as it is made of clear polyester and can be scratched by abrasive cleaners. Mild soap and water is recommended for general cleaning and disinfecting.

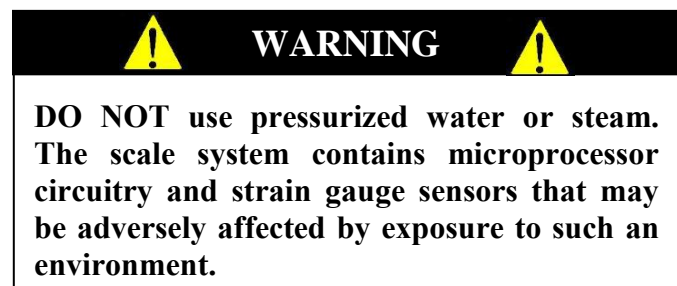

STORAGE

STORAGE

If storing this equipment for periods longer than three (3) months, remove the batteries. To maintain proper operation of this instrumentation, storage and transport conditions should not vary outside the following conditions: Relative Humidity 0% to 85%, Ambient Temperature 14°F to 122°F (-10°C to +50°C).

SPECIFICATIONS

MAXIMUM WEIGHT CAPACITY	1000 lb or 454 kg
PLATFORM SIZE	32 in x 40 in (81 cm x 102 cm)
DISPLAY TYPE	16-character Dot Matrix LCD
DISPLAY RESOLUTION	0.1 lb / 0.1 kg
ACCURACY	0.1% +/- 1 digit of displayed resolution for calibrated range
AUTO ZERO	One button operation
AUTO POWER DOWN	Approximately 30 seconds (adjustable to 300 seconds)
RECALL	Recalls last stored stable weight, height and BMI
AVERAGING	Automatic digital filter
POWER SUPPLY	Six (6) "D" cell batteries
CALIBRATION	Calibration is traceable to NIST standards
OPERATING CONDITIONS	Normal operating conditions for this product: Ambient Temperature Range: 68°F to 85°F (20°C to 30°C) Relative Humidity Range: 0% to 85% Avoid exposure to high-pressure water or steam.
TRANSPORT and STORAGE	Storage and transport conditions should not vary outside the following conditions: Relative Humidity 0% to 85%, Ambient Temperature 14°F to 122°F (-10°C to +50°C). Remove batteries if storing longer than three (3) months.

BUTTON FUNCTIONS



Figure 6: Button Display

ZERO / WEIGH



Press and hold to zero. The display will read “**WT = 0.0 Lb**”. This is used to zero the system before placing the patient on the scale system. This action also resets previously stored weight, height and BMI values to zero. Ensure that nothing is in contact with the weighing surface during this procedure.

Press to weigh. Weight stable indicator “**□**” flashes then remains solid when stable. Auto stores stable weight in memory.

SEND (PRINTER / EHR)



Press to send stored values to EHR or to printer. Output values include time, date, weight, height, and BMI. Also, indicated on display when paper is out, “**PAPER OUT**” and when door is open, “**PAPER DOOR OPEN**”.

RECALL



Press to recall the last stable weight. Display will scroll thru “**WT**”, “**HT**”, and “**BMI**” stored values in approximately two second intervals.

ENTER



Press to save height data and display BMI calculation. Stable weight, height and BMI values are stored in memory until next stable weight is stored or until system is zeroed.

BMI



Press to calculate the BMI. When the “**BMI**” is pressed, the default starting value “**HT = 65 in**” or “**HT = 165 cm**” is displayed. If there is no stored stable weight, the display will indicate “**NO WEIGHT DATA**” and then go back to the weigh screen “**WT = 0.0 Lb**”.

Continued next page

BUTTON FUNCTIONS cont'd

MENU



Press **Menu** to edit setup.

Setting the **UNITS**: Use **UP** or **DOWN** arrow buttons to select “**Lb**” or “**Kg**”. Press **ENTER** to save changes.

Setting **ON TIME**: Use **UP** or **DOWN** arrow buttons to adjust the “**ON TIME**”. The “**ON TIME**” may be set from 30 to 300 seconds in 30 second increments. Press **ENTER** to save changes.

Setting **TIME** and **DATE**: Use the **UP** arrow button to select digit. To change digit use the **DOWN** arrow button. Press **ENTER** to save changes.

NOTE: When selected, the year position defaults to “00”

UP



Press **UP** to adjust height up from the default, to increase the scale’s “on time”, or to select a digit when setting time and date.

DOWN



Press **DOWN** to adjust the height down from the default, to decrease the scale’s “on time”, or to change the value of a selected digit when setting time and date.

BASIC SYSTEM OPERATION

SETTING SYSTEM ZERO




Make sure the scale is free and clear of any obstructions and press and hold the **ZERO / WEIGH** button. The displayed message will indicate “**HOLD TO ZERO**” and count down to zero. Make sure that nothing is in contact with the scale while zeroing the system. In a few seconds, the display will read “**WT = 0.0 Lb**” (or **Kg**). This action also resets previously stored weight, height, and BMI values to zero.

**WARNING**
**CARE-GIVER MUST ASSIST WHEELCHAIR-BOUND PATIENTS.
DO NOT LEAVE PATIENT UNATTENDED ON THE SCALE PLATFORM.**

WEIGHING



Position the patient on the scale. The care-giver must be standing on the scale when moving the wheelchair on or off the scale. Lock the patient's wheelchair to prevent movement. Press the **ZERO/WEIGH** button. The weight stable indicator “” flashes on the display. When the weight is stable, the weight stable indicator remains solid. The display will indicate the patient's weight in either pounds or kilograms; example: “**WT = 123.5 Lb**”. The stable weight is auto stored in memory.

NOTE: If patient will be using a cane for support on the scale, place cane on scale while zeroing the system. This will ensure that the patient's **NET** weight will be displayed. It is recommended that the system be zeroed prior to each new patient.

CALCULATING BMI



Press to calculate BMI.

When the **BMI** button is pressed, the default starting value is set either to 65 in or to 165 cm; example: “**HT = 65 in**”. Use the **UP** or **DOWN** arrows to adjust the default height to the value of the patient's height.

NOTE: To calculate the BMI, a patient's stable weight needs be stored as indicated in “**WEIGHING**” above, if no stable weight has not been stored “**NO WEIGHT DATA**” will be displayed and the display will then go back to the weigh screen “**WT = 0.0 Lb**”.

Press **ENTER** to save height data and display BMI calculation. The BMI will be displayed; example: “**BMI = 20.5**”.

Continued next page

BASIC SYSTEM OPERATION con't


RECALLING LAST STABLE WEIGHT



Press to recall last stable weight.

The display will scroll thru “WT”, “HT” and “BMI” stored values in approximately two second intervals; example: “RCL WT = 123.5 Lb”, “RCL HT = 65 in” and “RCL BMI = 20.5”.

BATTERY REPLACEMENT

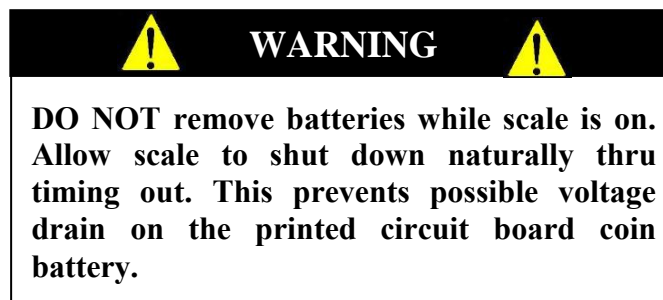

STEP 1: The display will read “REPLACE BATTERY”.

STEP 2: (**Figure 7**) Unscrew the panel screw on the Battery Compartment Cover and remove the battery compartment cover.

STEP 3: Remove and replace ALL six (6) “D” cell batteries.

STEP 4: Press the “WEIGH” button to confirm display is working.

STEP 5: Secure the battery cover using the panel screw.

STEP 6: Press and hold “WEIGH” button to zero the system.



Figure 7: Battery Replacement

THEORY OF OPERATION

SR Instruments patient weighing systems are digital scales. Strain-gauge force cells convert the force of an applied weight into an analog signal. This signal is amplified by an operational amplifier and converted to a digital signal by an on-chip analog to digital converter. The digital signal is filtered, converted to appropriate units, and displayed on the liquid crystal display.

Strain-gauge force cells each contain four strain gauges mounted in a full Wheatstone-bridge configuration. These bridges convert the physical movement of the force cell, due to the applied mass on the system, into minute changes in electrical resistance. These changes in resistance produce a voltage difference across the Wheatstone-bridge, which is amplified by the operational amplifier. The amplifier is configured to current sum the output of each cell, with potentiometers serving to normalize the sensitivity (voltage out per unit of weight applied) of each bridge. The offset potentiometer produces a small current, which nulls the output of the amplifier for an unloaded system.

The output of the operational amplifier is digitized by the analog to digital converter. The sigma-delta converter sums a rapid sequence of 0's (0 volts) and 1's (reference voltage) to achieve balance with the input signal from the amplifier.

The micro-controller filters the digital output of the analog to digital converter, subtracts the value saved during the system zero operation and scales the filtered output, and then displays the result on the liquid crystal display. The micro-controller performs a moving-median filter of data for continuous weigh the micro-controller performs checks for signal stability before locking in on the reading.

The micro-controller can be placed in a calibration mode, where the system can be re-calibrated. In the calibration mode, the system slope is calculated from two points (zero and full scale) in the 2-point calibration mode or the slope and change in slope is calculated from three points (zero, half, and full scale) in the 3-point calibration mode.

CALIBRATION

NOTE: Ensure that nothing is in contact with the scale system during this procedure. Remove hands from the system when noting the displayed calibration results.

CHECKING CALIBRATION

STEP 1: Select known calibrated weights, traceable to NIST.

DO NOT USE barbells or uncalibrated weights.

STEP 2: Zero the scale by pressing and holding the **ZERO / WEIGH** button.

STEP 3: Place the calibrated weight on the scale. Wait for scale to stabilize; note scale reading.

STEP 4: Scale readings for both weights should be within Calibration Tolerance Table (**Figure 8**)

CALIBRATION TOLERANCE TABLE		
LOW LIMIT	APPLIED LOAD	HIGH LIMIT
99.9	100.0	100.1
199.8	200.0	200.2
299.7	300.0	300.3
399.6	400.0	400.4
499.5	500.0	500.5
599.4	600.0	600.6
999.0	1000.0	1001.0

Figure 8: Calibration Tolerance Table

⚠ IMPORTANT ⚠

CALIBRATION Qualified service personnel only should perform this procedure. The SR735i load cells have no user serviceable components and should not be tampered with for any reason. Re-calibration is generally not required, but should be verified periodically to ensure accuracy. The recommendation for calibration check is at least once every 12 months, or as individual maintenance policy requires.

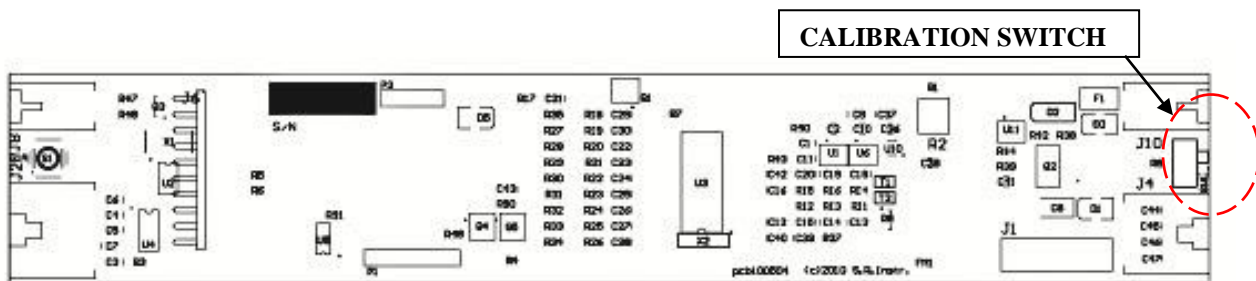


Figure 9: Calibration Switch Diagram

⚠ CAUTION ⚠

The integrated circuits and semiconductors on the printed circuit boards may be damaged by electrostatic discharge (ESD). Be sure to use proper handling precautions at all times.

Continued next page

CALIBRATION cont'd

SETTING CALIBRATION

NOTE: Ensure that nothing is in contact with the scale system during this procedure. Remove hands from the system when noting the displayed calibration results.

STEP 1: Remove the three (3) screws on the left end cover of the display housing. Remove the cover. Put the scale system into the Calibration Mode by switching the calibration switch on the display board (**Figure 9**).

STEP 2: Select known calibrated weights, traceable to NIST, up to the Full Scale value (maximum capacity).

STEP 3: Press the “MENU” button until “1000Lb FULL” is displayed. Set the FULL value to the actual quantity of calibrated weight being used for Full Scale. Use the “UP” arrow button to select the digit and the “DOWN” arrow button to change digit. Press “ENTER” to save changes. Press “WEIGH” button to abort any changes.

STEP 4: Press the MENU button until “500Lb HALF” is displayed.

NOTE: The Half-Scale value is a value between zero and the Full Scale values. It is usually close to half the Full Scale value.

Set the HALF value to the actual quantity of calibrated weight being used for Half Scale. Use the “UP” arrow button to select the digit and the “DOWN” arrow button to change the digit. Press “ENTER” to save changes. Press “WEIGH” to abort the changes.

STEP 5: Press MENU button until “3 PT CAL” is displayed. Press the “UP” arrow button.

STEP 6: Zero the scale by removing all weight from the platform. Press the “UP” arrow button to “ADD HALF”.

STEP 7: Place the HALF calibrated weight on platform. Allow weight to stabilize. Press “UP” arrow button to save change.

STEP 8: Place the FULL calibrated weight on platform. Allow weight to stabilize. Press “UP” arrow to save change.

STEP 9: Press “ENTER” button to save the calibration, or “WEIGH” button to exit without saving.

STEP 10: Switch the scale system out of the Calibration Mode on the display board (**Figure 9**).

STEP 11: Place the left end cover on the display housing. Attach using three (3) cover screws.

CONFIGURING SEND BUTTON

STEP 1: Remove the three (3) screws on the left end cover of the display housing. Remove the cover. Put the scale system into the Calibration Mode by switching the calibration switch on the display board (**Figure 9**).

STEP 2: Press the “MENU” key until “DATA OUT =PRINT” is displayed

STEP 3: Press “UP” or “DOWN” button to select between “DATA OUT =NONE”, “DATA OUT =PRINT” or “DATA OUT =EHR”

STEP 4: Press the “ENTER” button to save.

STEP 5: Switch the scale system out of the Calibration Mode on the display board (**Figure 9**).

STEP 6: Place the left end cover on the display housing. Attach using three (3) cover screws.

EHR OUTPUT

Electronic Medical Records (EMR) / Electronic Health Records (EHR) technology is supported by the i-series scales from SR Instruments. The EMR/ EHR software is not included with scale. There are many different EMR/ EHR software available and the connectivity of the scale to the software requires the services of a professional.

A null modem cable is required to communicate with the scale.

Scale output when send button is configured to Electronic Health Records is as follows:

```
<esc>R<esc>Wnnn.n<esc>Uuu<esc>E
```

Where: R is read

W is weight

nnn.n is the weight in Lb or kg

U is units (lb or kg; in or cm)

uu is LB or KG

E is end of packet.

TROUBLESHOOTING

SYMPTOM	REASON/CORRECTIVE ACTION
The characters only appear on half of the display.	Press the “ WEIGH ” button or remove one battery. Wait five seconds, then re-install the battery and try the “ WEIGH ” button again.
The display lights appear to work, but do not respond to button activation.	Check to ensure the membrane switch label is correctly plugged into display board. Check to ensure the calibration switch is not in the Calibration Mode (Figure 9).
The display shows no reading at all.	Check to ensure batteries are installed correctly (see directions for BATTERY REPLACEMENT). Check display cable to make sure it is connected securely.
For additional information or assistance, telephone our Service Hotline: 1-800-654-6360 or e-mail: sri@srinstruments.com	

WARRANTY

FOUR YEAR LIMITED WARRANTY

SR Instruments, Inc. systems are manufactured with high quality components. SR Instruments, Inc. warrants that all new equipment will be free from defects in material or workmanship, under normal use and service, for a period of four (4) years from the date of purchase by the original purchaser. Normal wear and tear, injury by natural forces, user neglect, and purposeful destruction are not covered by this warranty. Warranty service must be performed by the factory or an authorized repair station. Service provided on equipment returned to the factory or authorized repair station includes labor to replace defective parts. Goods returned must be shipped with transportation and/or broker charges prepaid. SR Instruments, Inc.'s obligation is limited to replacement of parts that have been so returned and are disclosed to SR Instruments, Inc.'s satisfaction to be defective. The provisions of this warranty clause are in lieu of all other warranties, expressed or implied, and of all other obligations or liabilities on SR Instruments, Inc.'s part, and it neither assumes nor authorizes any other person to assume for SR Instruments, Inc. any other liabilities in connection with the sale of said articles. In no event shall SR Instruments, Inc. be liable for any subsequent or special damages. Any misuse, improper installation, or tampering, shall void this warranty.

DAMAGED SHIPMENTS

Title passes to purchaser upon delivery to Transportation Company. Any claims for shortage or damage should be filed with the delivery carrier by purchaser.

RETURN POLICY

All products being returned to SR Instruments, Inc. require a Return Goods Authorization number (RGA). To receive an RGA, call our Technical Service Team at 716-693-5977 or toll-free in the USA and Canada at 800-654-6360.

When inquiry is made, please supply model and serial numbers, purchase order, if the scale was bought on contract, and reason for return.

Generally, deleted, damaged, and outdated merchandise will not be accepted for credit. A minimum restocking charge of 15% will be assessed on return of current merchandise.

All returns are to be shipped **FREIGHT PREPAID** to: SR Instruments, Inc., 600 Young Street, Tonawanda, NY 14150.

RESTOCKING FEE

- **15% fee** for any scale that has been opened and used
- **10% fee** for any scale returned that has been ordered incorrectly or refused delivery with no model change
- **5% fee** if an error in ordering has been made and a different model exchanged
- **No fees** will be charged if the scale is returned because of an error on the part of SR Instruments, Inc.
- **No returns** accepted after 60 days.

SR Scales[®]

by **SR**[®]Instruments, Inc.

**Precision & Technology in
Perfect Balance**[®]