

Quality. Reliability. Accuracy.

Visiting Nurse Scale Model SR241

The SR241 Portable Infant/Adult Scale is a unique portable system that offers the versatility of weighing infants, toddlers, and adults on one lightweight and moistureresistant system. Weighing in at a mere seven pounds, this scale is perfect for visiting nurse applications. Its sturdy construction and inherent durability also lends itself well for use in pediatrician or physician offices.



LIGHTWEIGHT

Weighs only 7 pounds. Built-in carrying handle.

HEAVY-DUTY

Injection molded out of high strength, impact-resistant plastic.

WEATHER TIGHT

The entire scale is designed with the visiting nurse in mind. So a little rain or snow is just fine.

EASILY CLEANED

There are no hidden crevices or dirtcatching seams. The entire scale can easily be cleaned to mitigate crosscontamination issues and ensure caregiver safety as well.

2-4-1

2 scales 4 the price of 1. Weigh both adults and infants on the same scale. Weigh babies up to 44 pounds and adults up to 400 pounds. The removable tray separates the adult-weighing surface from the baby-weighing surface.

BATTERY POWER

A single 9vDC provides up to 3000 weigh cycles.

OPTIONAL EQUIPMENT

Carrying case

Adult Mode



WARRANTY

THE USA

SR[®]Instruments.inc.

Platform size: 14" x 12"

600 Young Street, Tonawanda, New York 14150 USA T: 716.693.5977 · F: 716.693.5854 · TF: 1.800.654.6360 www.srscales.com email: info-request@srinstruments.com

Visiting Nurse Scale Model SR241





SYSTEM **SPECIFICATIONS MAX WEIGHT** 400 Lb /181 Kg **DISPLAY TYPE** 4-1/2 digit 7 segment LCD Adult Mode: 0.1 Lb (0.1 Kg) **DISPLAY RESOLUTION** Infant Mode: 1/4 oz (0.01 Kg) ACCURACY 0.2% + - 1 digit of displayed resolution **DISPLAY UNITS** Pounds / Ounces and Kilograms 1 common 9v Alkaline Battery. **POWER** Provides approximately 3000 weigh cycles. Approximately 60 sec. to conserve **AUTO POWER DOWN** battery life. A low battery indication will appear when it's time to replace (APD) the batteries. 4 Lb (1.8 Kg) (scale only) WEIGHT 7 Lb (3.2 Kg) (scale with cradle) WARRANTY 4 years **ACCESSORIES AND OPTIONS** DESCRIPTION PART NUMBER **INFANT / ADULT SCALE** SR241 **KILOGRAMS ONLY** SR241KG ADULT SCALE ONLY (NO TRAY) SR241S0 **KILOGRAMS ONLY** SR241SOKG ANTIMICROBIAL INFUSED SRC-600 PEDIATRIC SCALE CART **REPLACEMENT INFANT TRAY** SM4028 SOFT SIDED CARRYING CASE PP4022






DIMENSIONS in [cm]

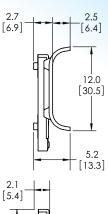


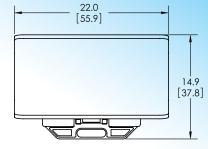


SCR

/ SR Instruments, Inc

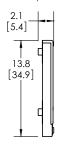

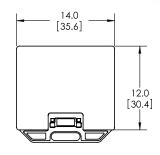

PACKAGING

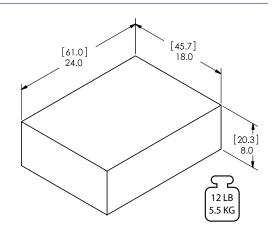

- Weather tight design
- No place for dirt to hide
- Easy to clean