

# Stand-On Scale

## Model SR585i

Now featuring touchless weighing\*

The SR585i Stand-On Scale is equipped with a large platform and full wrap-around handrail for patient stability.

Constructed of heavy-duty cast aluminum, the SR585i has a capacity of 1000 Lbs, yet is lightweight and easily transported by caregivers.

Image shown with optional handrails and height bar



Platform size: 18" x 20"

\*SR555i, SR585i, SR755i scales now feature optional AC power that enables scales to remain plugged in, always on, and ready to weigh a patient with no button activation.



### ANTIMICROBIAL

SR Scales exclusive: contains Microban<sup>®</sup> antimicrobial technology for powerful, continuous product protection against bacteria.

### PORTABLE

Just Tip n' Go! Take the scale to the patient for an ideal, point-of-use scale.

### MEMORY RECALL

Recalls the last stored weight.

### AC OR BATTERY OPERATED

An optional medical grade wall power supply can be plugged into the scale for touchless scale operation. Additionally 6 "D" cell batteries provide approximately 10,000 weigh cycles with programmable Auto Power Down (APD) to conserve power when not plugged in.

### SOLID CAST ALUMINUM CONSTRUCTION

Durable and easy-to-clean powder coated finish.

### WIDE GRIP, WRAP-AROUND HANDRAILS

Patient can use the wrap-around handrails and integrated front-hand hold for stability without affecting accuracy.

### BMI CALCULATION

Enter height information to calculate and display Body Mass Index.

### EHR/EMR READY

Serial data output can be used to communicate with electronic health records through send or get commands.

### OPTIONAL EQUIPMENT

- Medical grade wall power supply
- Thermal printer
- Height Rod - Easily measure in both inches and centimeters.



EASY TRANSPORT  
& STORAGE



4 YEAR  
WARRANTY

# Stand-On Scale

## Model SR585i



### SYSTEM SPECIFICATIONS

|                              |  |
|------------------------------|--|
| <b>MAX WEIGHT</b>            | 1000 Lb / 454 Kg   |
| <b>DISPLAY TYPE</b>          | 16 Character LCD w/ LED Backlight  |
| <b>DISPLAY RESOLUTION</b>    | 0.1 Lb/Kg  |
| <b>ACCURACY</b>              | 0.1% +/- 1 digit of displayed resolution   |
| <b>DISPLAY UNITS</b>         | Pounds / Kilograms   |
| <b>POWER</b>                 | 6 "D" cell alkaline batteries and optional medical grade wall power supply   |
| <b>BATTERY LIFE</b>          | 10,000 Weigh Cycles (est. based on 30 second on-time)  |
| <b>AUTO POWER DOWN (APD)</b> | Programmable on-time to conserve battery life. A low battery indication will appear when it's time to replace the batteries. |
| <b>MEMORY</b>                | Last stable weight is automatically stored in memory. Press RECALL to retrieve.  |

### BMI CALCULATION

Press BMI and input height. Press ENTER to view BMI information.

### SERIAL DATA OUTPUT (EHR)



9-pin serial data output can be used to send weight data to optional printer or EHR/EMR Systems.

### ANTIMICROBIAL


See [www.microban.com/sr-scales](http://www.microban.com/sr-scales) for specifications on Microban® technology.

**WARRANTY** 4 years


### ACCESSORIES AND OPTIONS

| DESCRIPTION   | PART NUMBER   |
|---|---------------|
| <b>SR585I WITH AC ADAPTOR</b>   | SR585i-AC     |
| <b>KILOGRAM ONLY SCALE</b>  | SR585i-KG     |
| <b>SR585I WITH HEIGHT ROD</b>   | SR585i-H      |
| <b>PRINTER KIT</b><br>includes FRAP1300, FRTP130012C, FRBC1300, SE8294, and all mounting hardware | SR8220        |
| <b>REPLACEMENT HEIGHT ROD</b><br>(ADD CM SUFFIX FOR CENTIMETER)                                   | SR9110        |
| <b>MEDICAL GRADE POWER SUPPLY</b>   | FRGSM06U09P1J |

### DIMENSIONS in [cm]

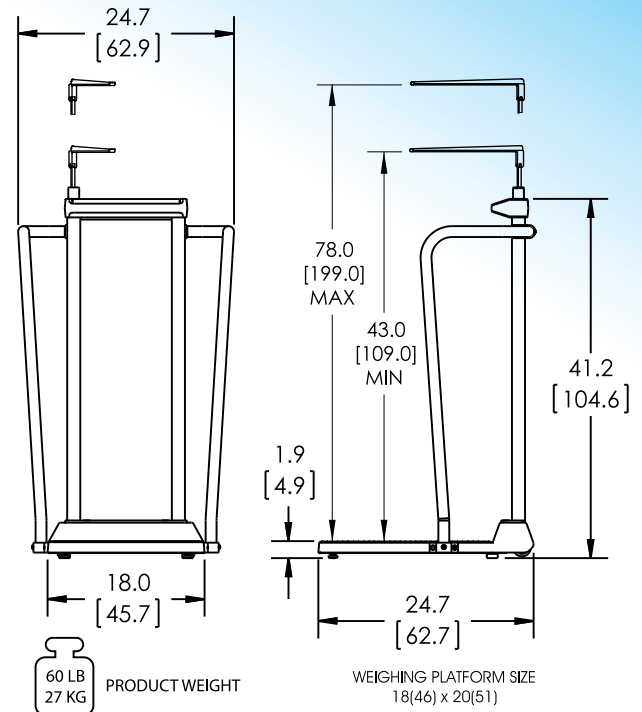

### PACKAGING

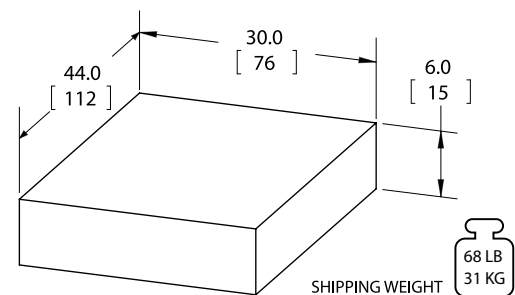

**Integrated front hand hold and wrap-around handles promote patient stability.**