

SR Scales®

"The Measure of Excellence"

Heavy Duty Scale Model **SR** V942

When you need a scale that is durable and heavy duty, consider the SRV942. While used primarily in stationary veterinary practices, this scale also has a wide range of application in the zoological field.

Constructed of stainless steel with a weighing capacity of 550 pounds, the SRV942 is extremely durable and provides consistent weight readings.

The rechargeable battery permits remote weighing when required. The SRV942 is equipped with a wall mount and power adaptor for stationary applications.

This scale is equipped with wheels and a fixed handle to allow for easy transport when required.

Whatever your weighing needs may be, **SR Scales®** is *"The measure of Excellence."*

Stainless Steel Platform

Durable and Easy to Clean

Surface Mat

Removable for Quick and Easy Cleaning

Platform Size- 22"X 36"

Auto-Power Shutdown

Automatically Shuts down after Two Minutes

Multiple Power Source Options

Rechargeable battery with AC/DC Adaptor



**2
YEAR
WARRANTY**

SR® Instruments, Inc. 600 Young Street, Tonawanda, NY 14150 USA
TF: 1.888.WEIGH.ME Ph: 716.693.5977 F: 716.693.5854
www.srscales.com email: info@srscales.com

**SR®
inside**


SYSTEM SPECIFICATIONS

Max Weight	550lb or 249kg
Display Resolution	.2 lb or .1kg
Display Units	Pounds/ Kilograms
Four Buttons	On/Zero/Off, Tare, Unit and Hold
Hold	Freezes the weight on the display
Tare	Allows objects to be placed on scale and zeroed out
Power	Rechargeable battery with DC Adaptor
Surface Mat	Removable for easy cleaning
Portability	Built-in handle and wheels
Auto- Power Down	Approximately two (2) minutes of non-use
Warranty	Two (2) Years

Heavy Duty Scale

Model **SR** V942

DIMENSIONS in (cm)

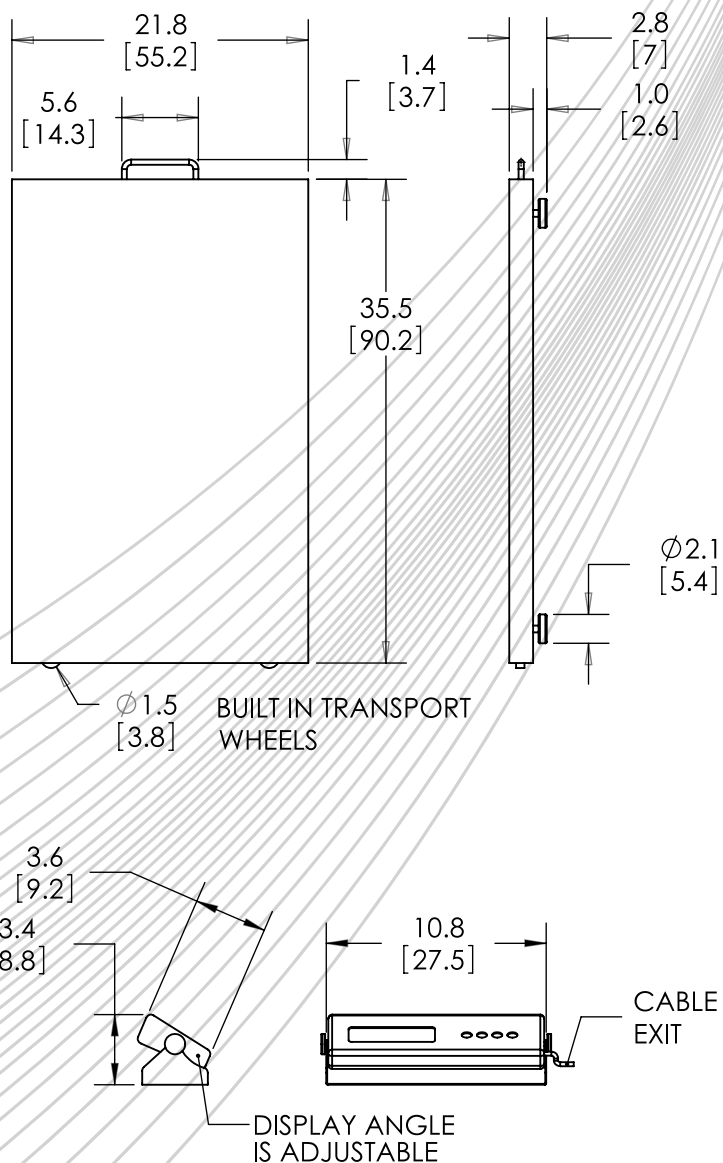

ACCESSORIES and OPTIONS

Description	Part Number
Replacement Mat	SRV942M
Replacement Charger	SRV942CH

PACKAGING

